



Dimensions

Unit : mm
Tolerance in dimensions : ±10%


■ Unit Mount

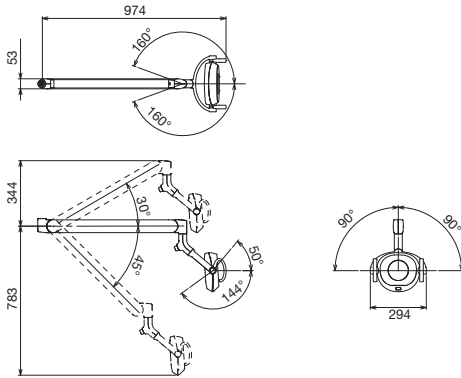

■ Post Mount

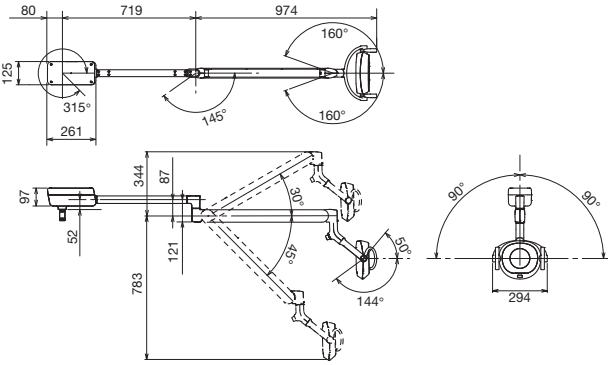


■ Ceiling Mount

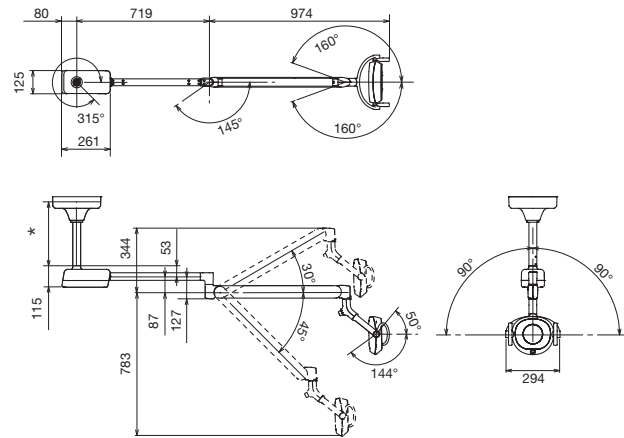

*Select suspension tube length based upon ceiling height.

■ Track Mount

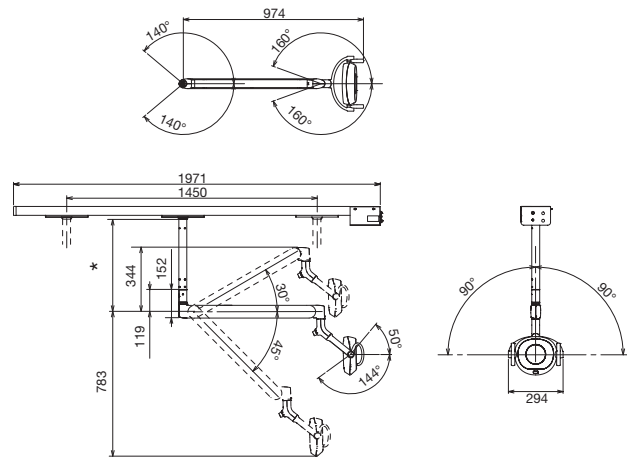

*Select suspension tube length based upon ceiling height.

EURUS LIGHT



Brightly & Gently

EURUS LIGHT—the next-generation LED lighting fixture (conforming to ISO9680), capable of giving off light with variable intensity from 5,000 to 34,000 luxes, reproduces natural light with a high functional color rendering capacity. With the touchless sensor, the light can be turned on and off and the composite mode can be adjusted.



Installation Variations



Unit Mount



Post Mount



Ceiling Mount



Track Mount



Flexible lighting

The Belmont spring arm allows the light to be moved around and stay exactly where it is desired so that the dentist/dental assistant can start treatment without delay. The arm covers a wide area, including the maxillary and mandibular areas of the patient, for sufficient illumination.



Stay clean

Lens cover and autoclavable light handle can be removed with a single lever action for easy cleaning.

Specifications

Rated voltage	DC 19 V
Rated input.....	1.3 A
Operation mode	Continuous operation
Light source	6 LED lamps
Optical performance	Standard irradiation distance: 650 mm
	■ In treatment mode Central illuminance: 5,000–34,000 lx Correlated color temperature: 5,000 K
	■ In Composite Safe mode (Compliant with ISO 9680: 2014 5.2.10) Central illuminance: 11,000 lx Correlated color temperature: 2,700 K
Illuminance pattern	Pattern dimensions (length x width): 85 x 155 mm
Color rendering index	More than 90

Options



Patient mirror